


# EEGLAB STEPS

---

Hayat Chamtie




# 1) Call EEGLAB in MATLAB

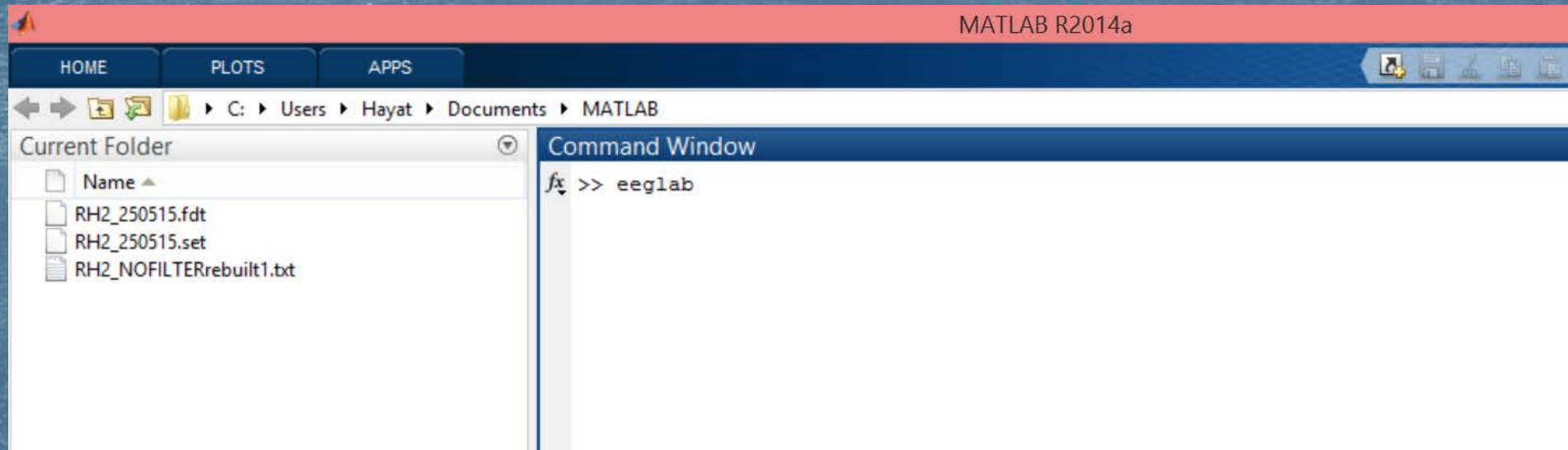

Note: Make sure to place EEGLAB file in Matlab folder and set the path to it in order for it to be called.




# Import Biosemi file

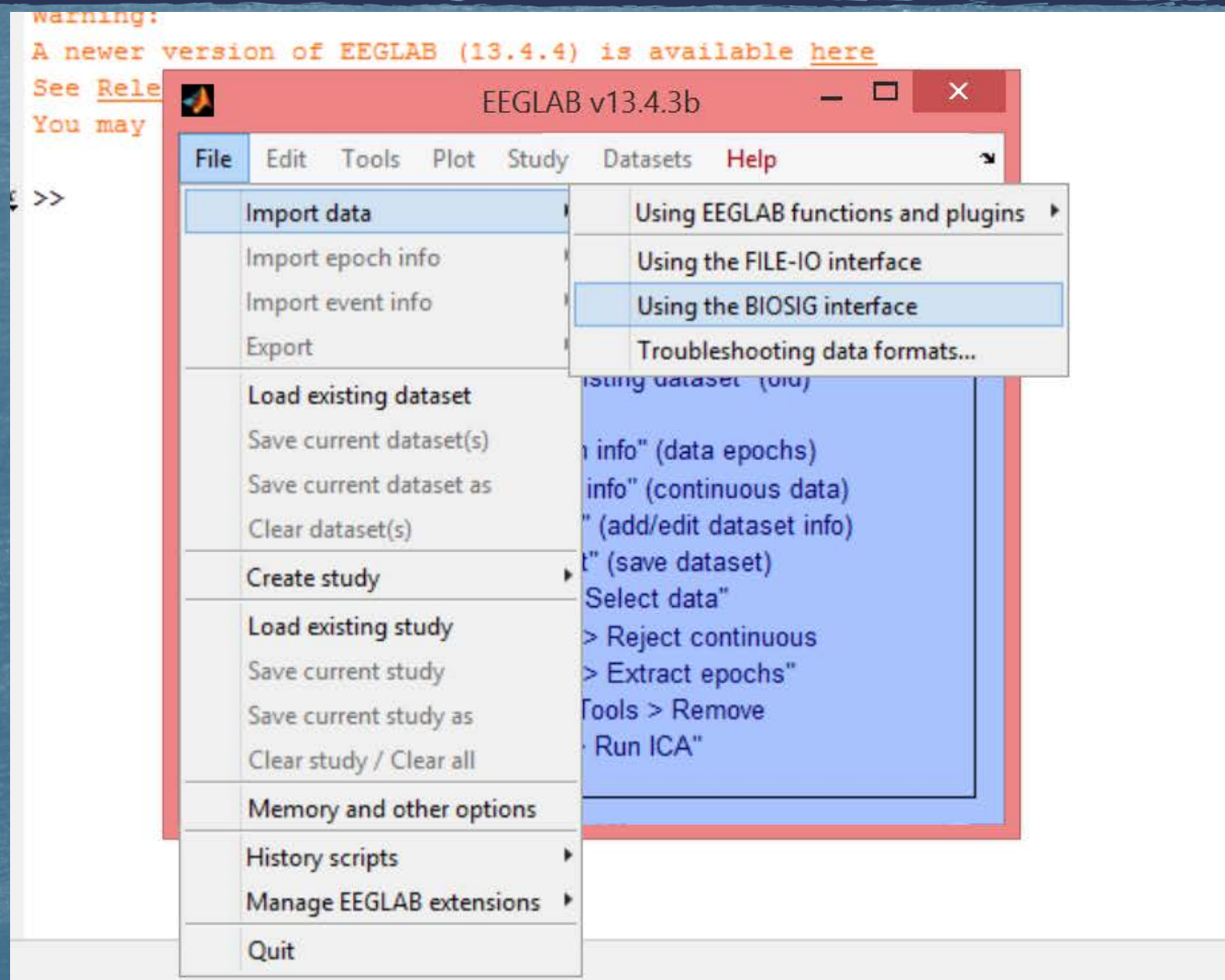




Change sampling rate to 512Hz

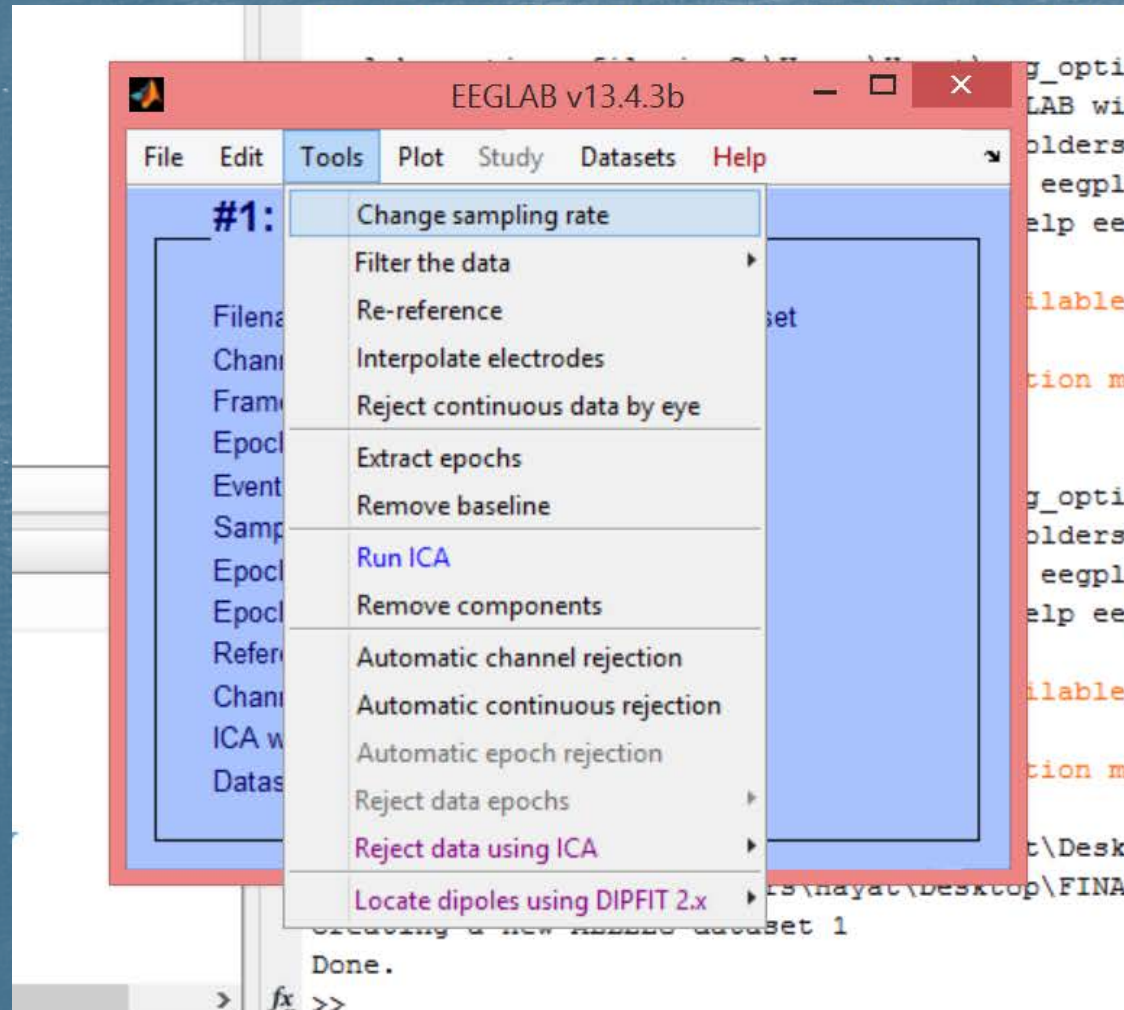

Filter the data  
Do each of the three as a separate step



- 1) Low pass : 1
- 2) High pass: 60
- 3) Notch filter: Low-49, high-51

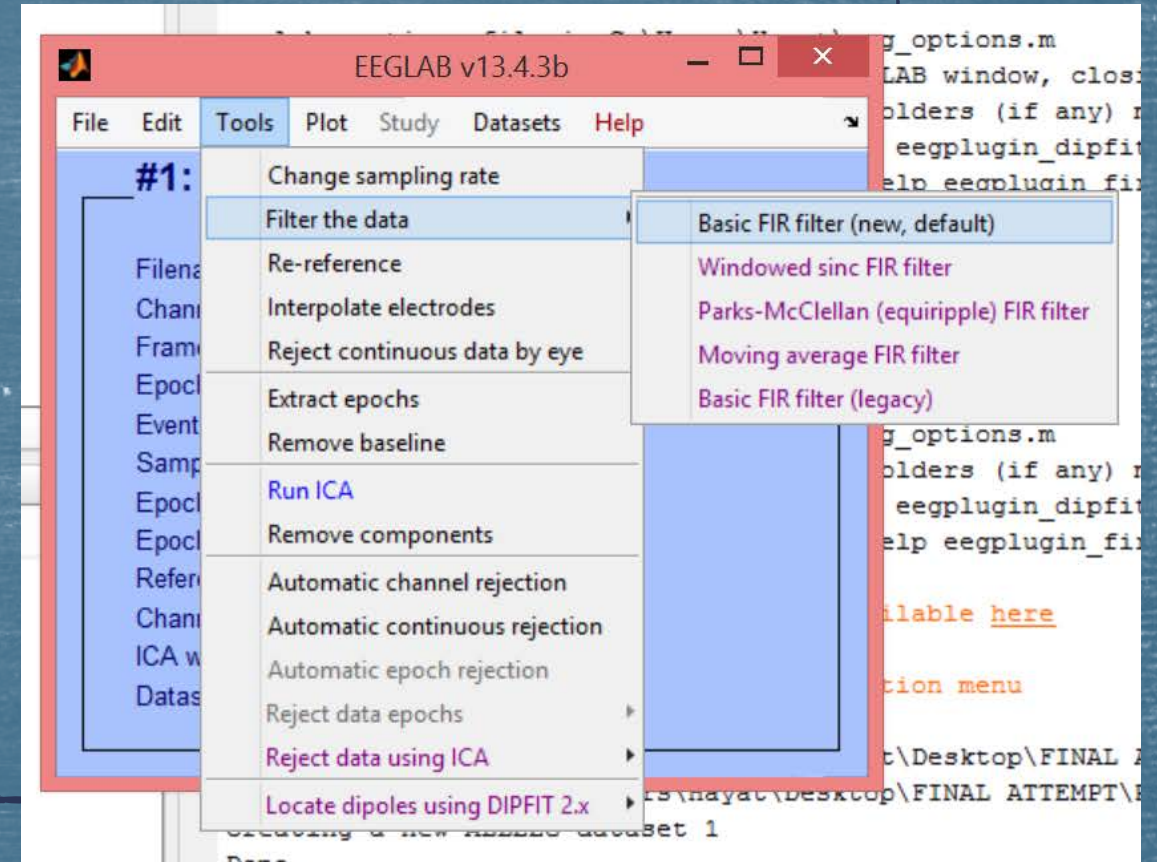




# Deleting extra channels

Tick the box below and remove all channels after EXG3 (including it)

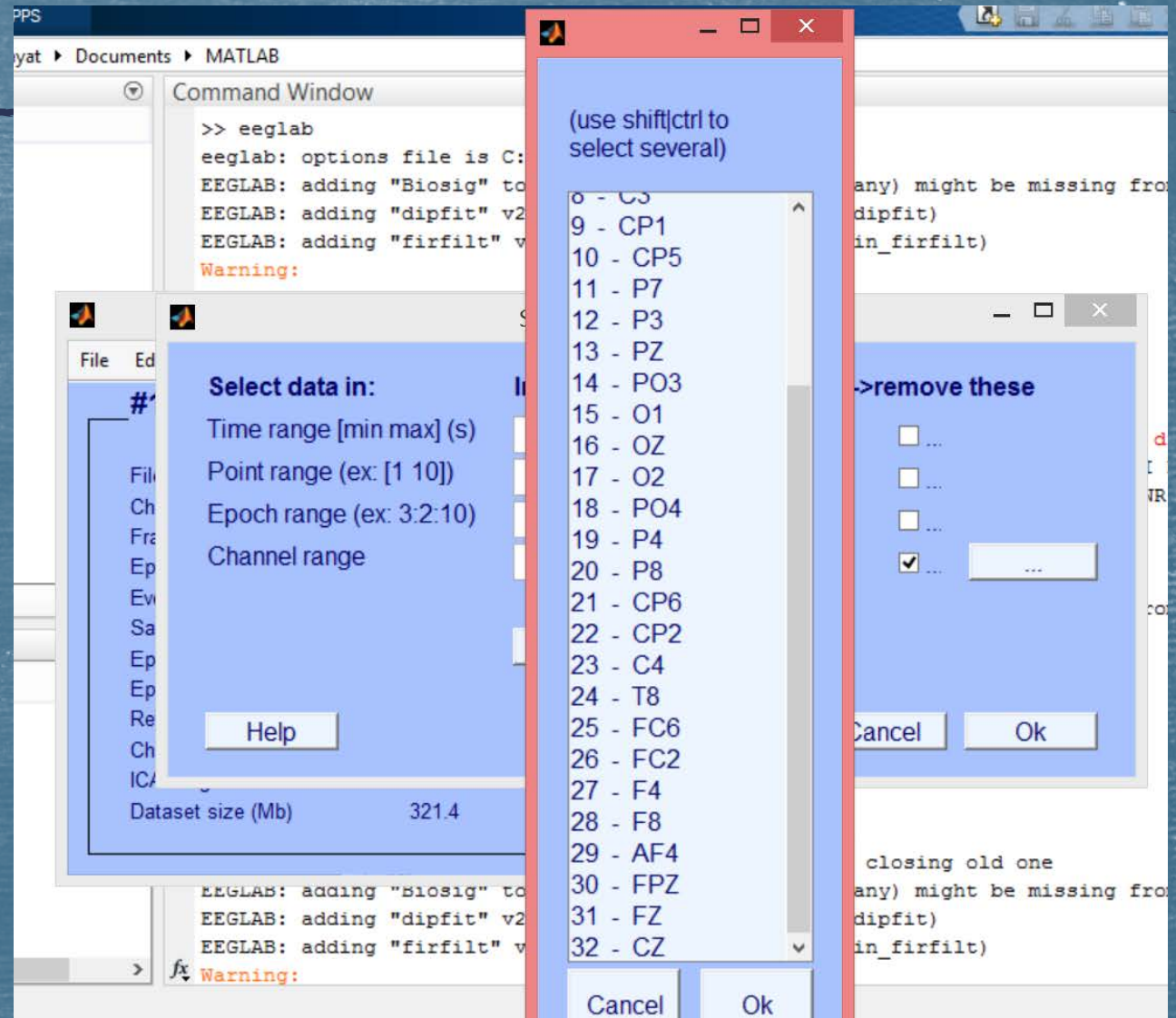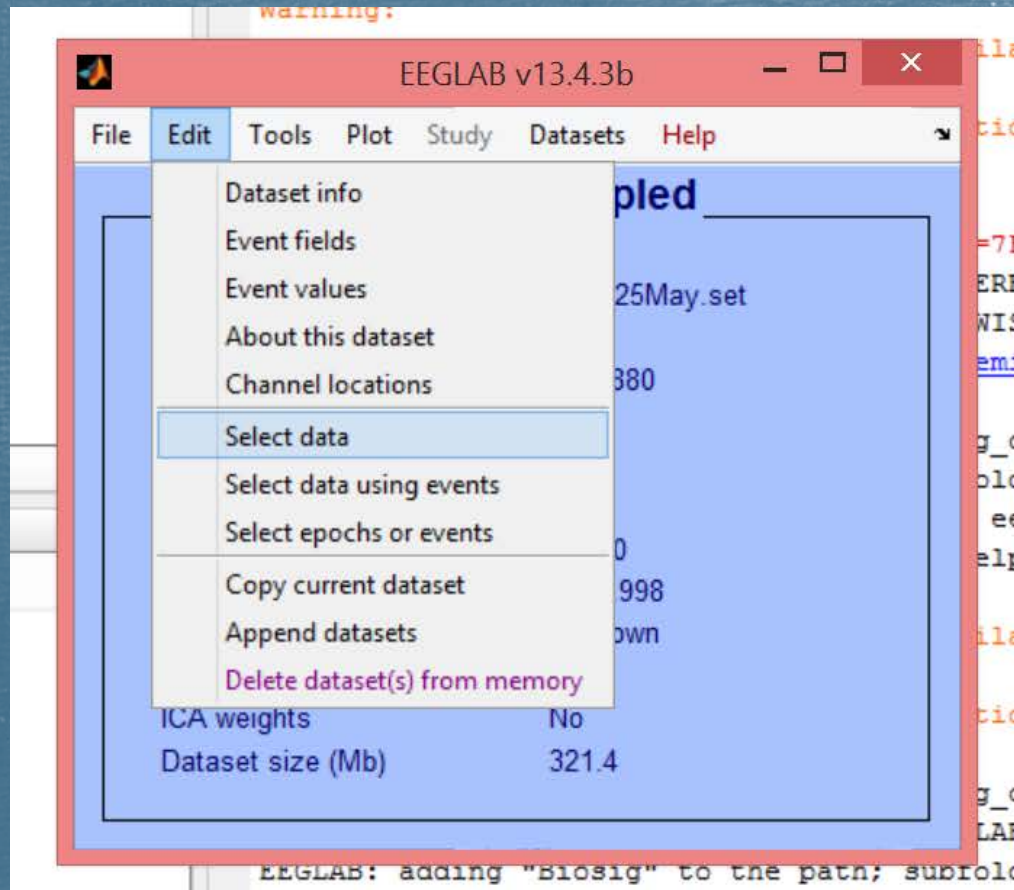




# Referencing of EXG1 and EXG2

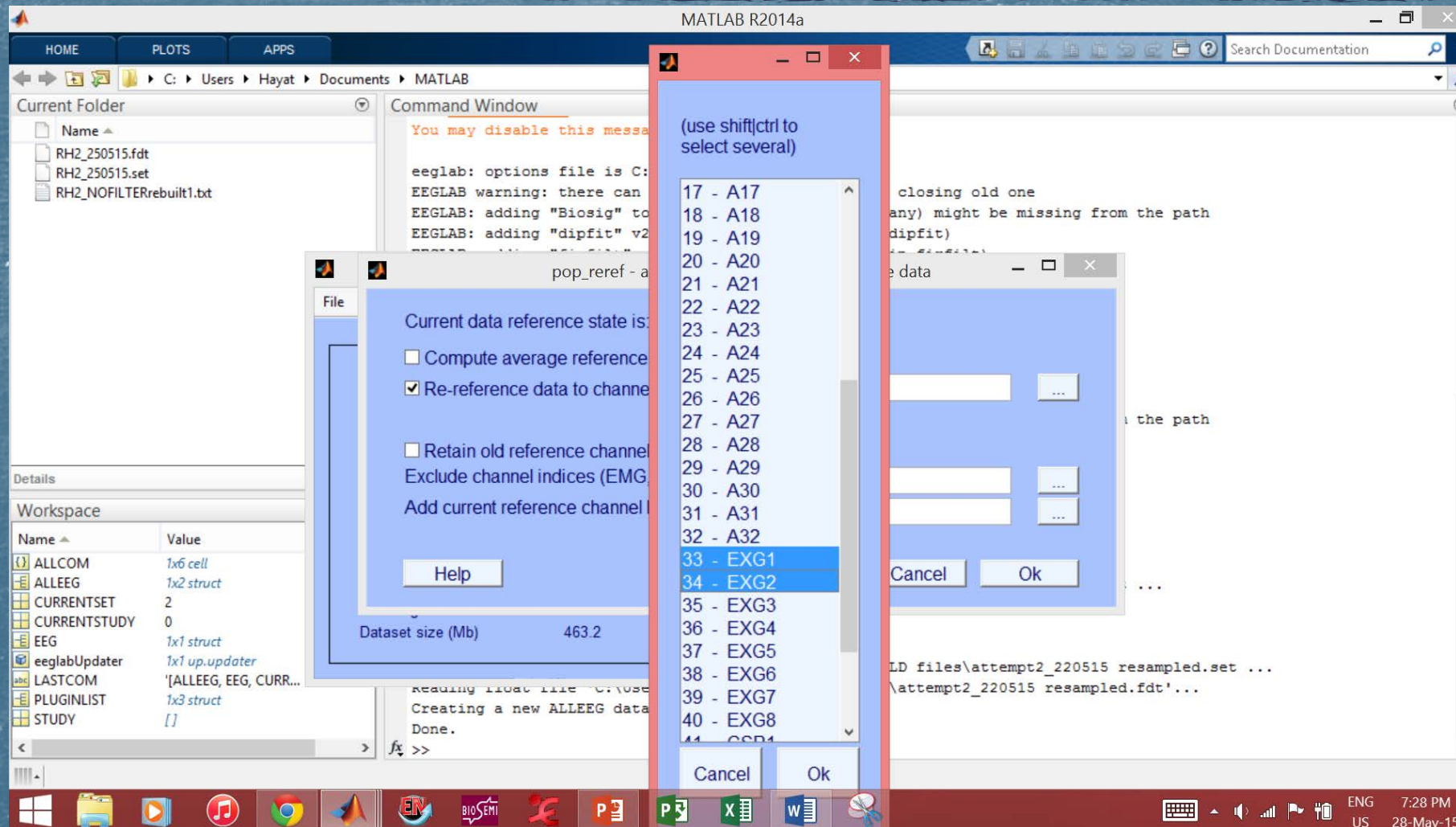




# Name channels

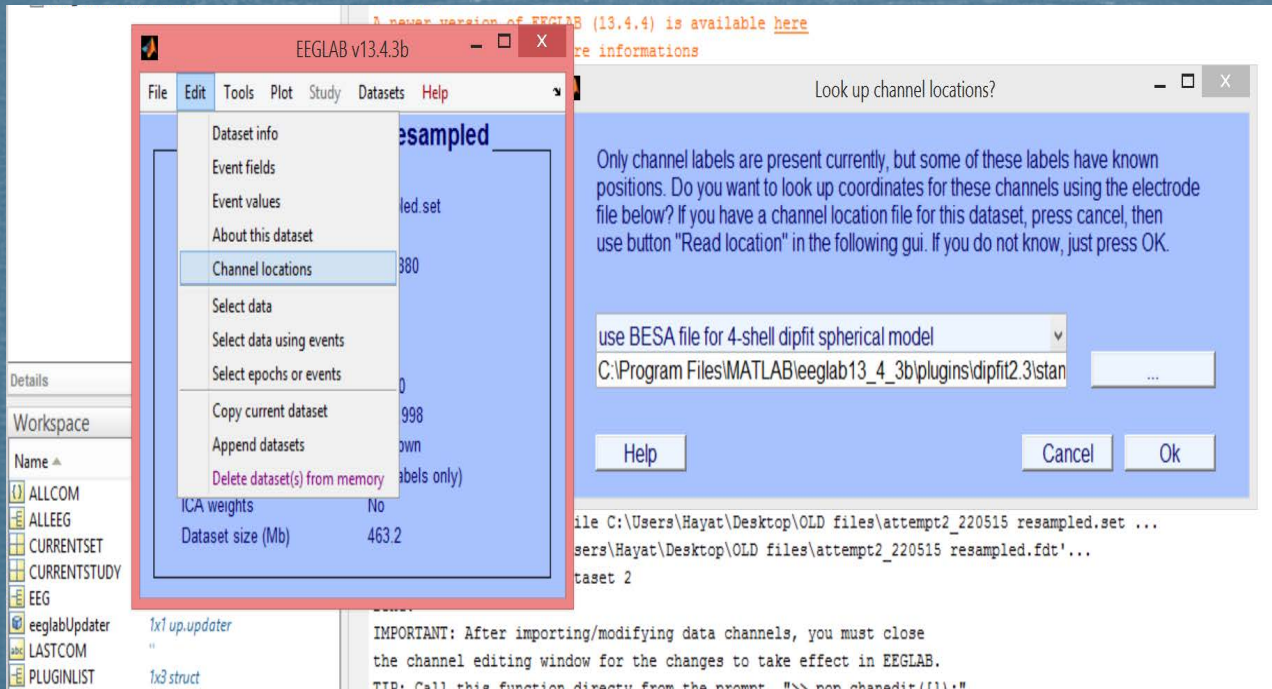


- ▶ Click 'cancel'
- ▶ Read channel locations from 'Biosemi32.ced' file

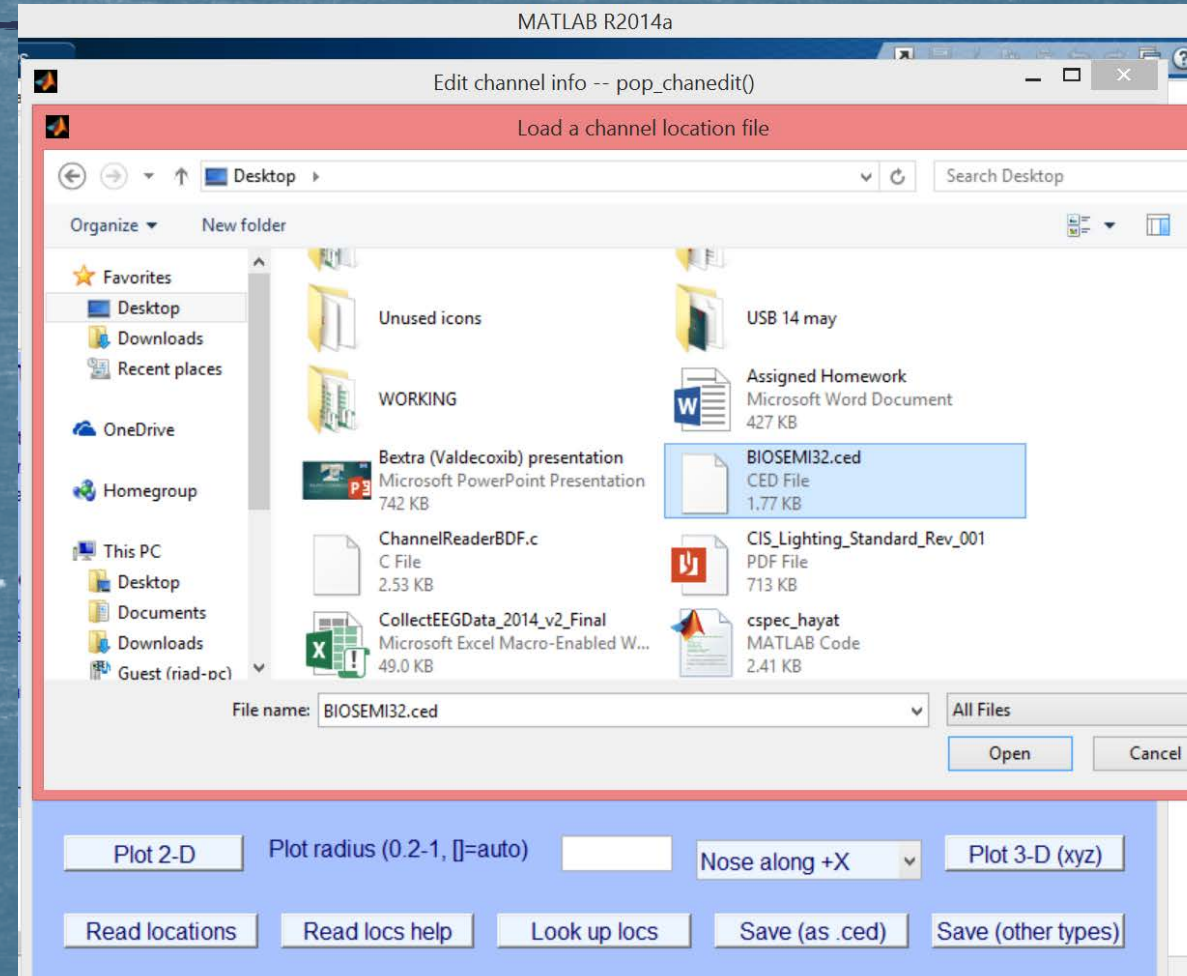




# Exporting to text file

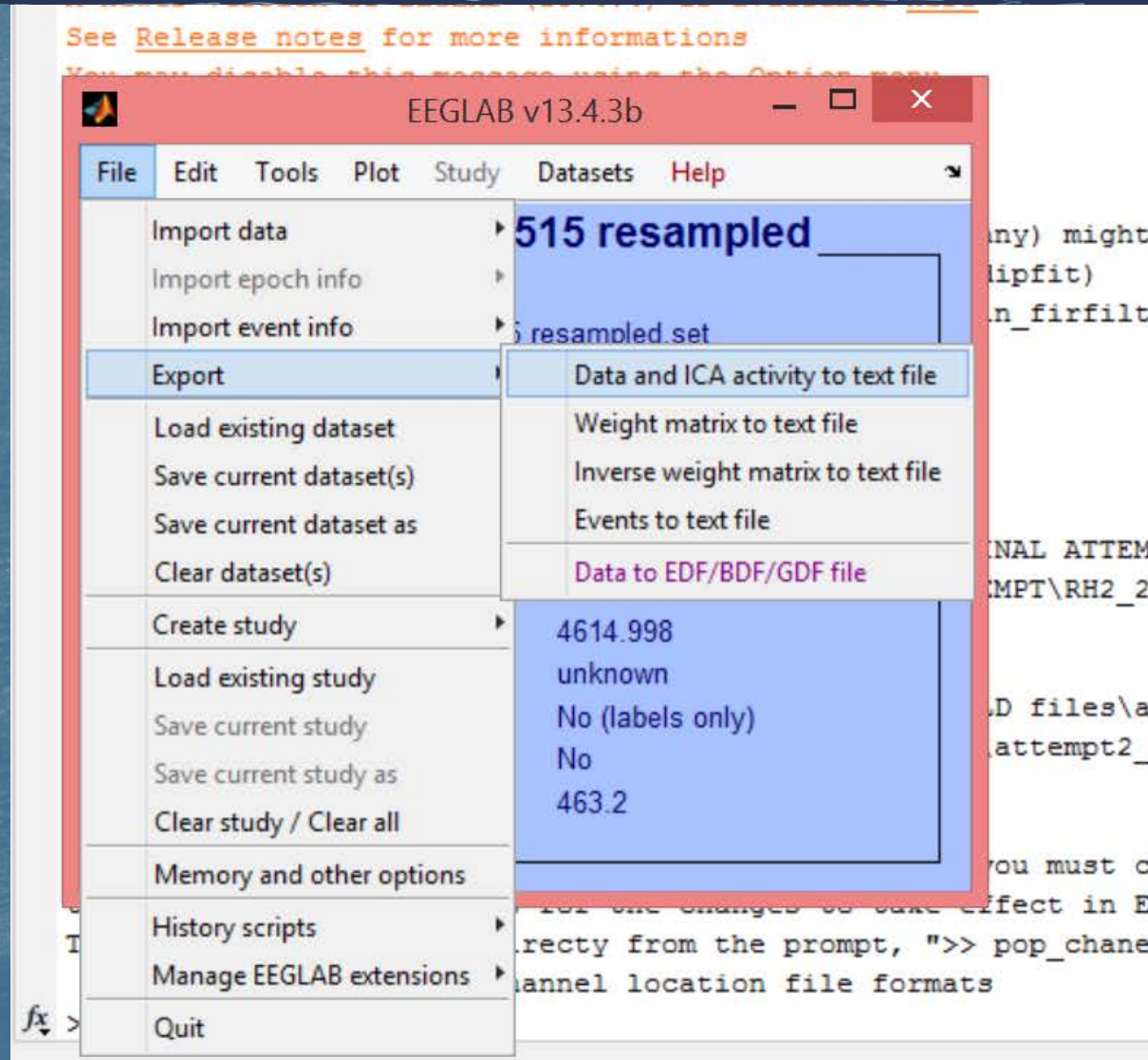





# Performing ICA analysis & removing Eye blinks, ECG and EMG artifacts

To be done after data is cut into smaller files before power spectra and coherence correlation analysis

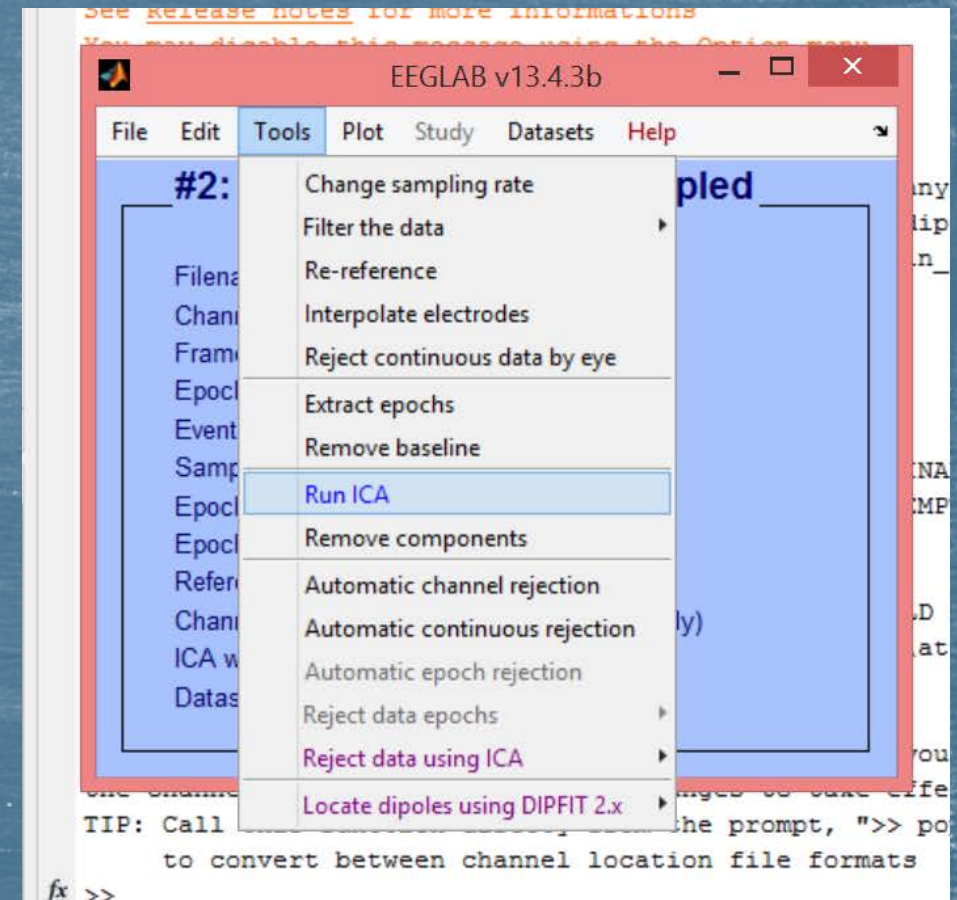



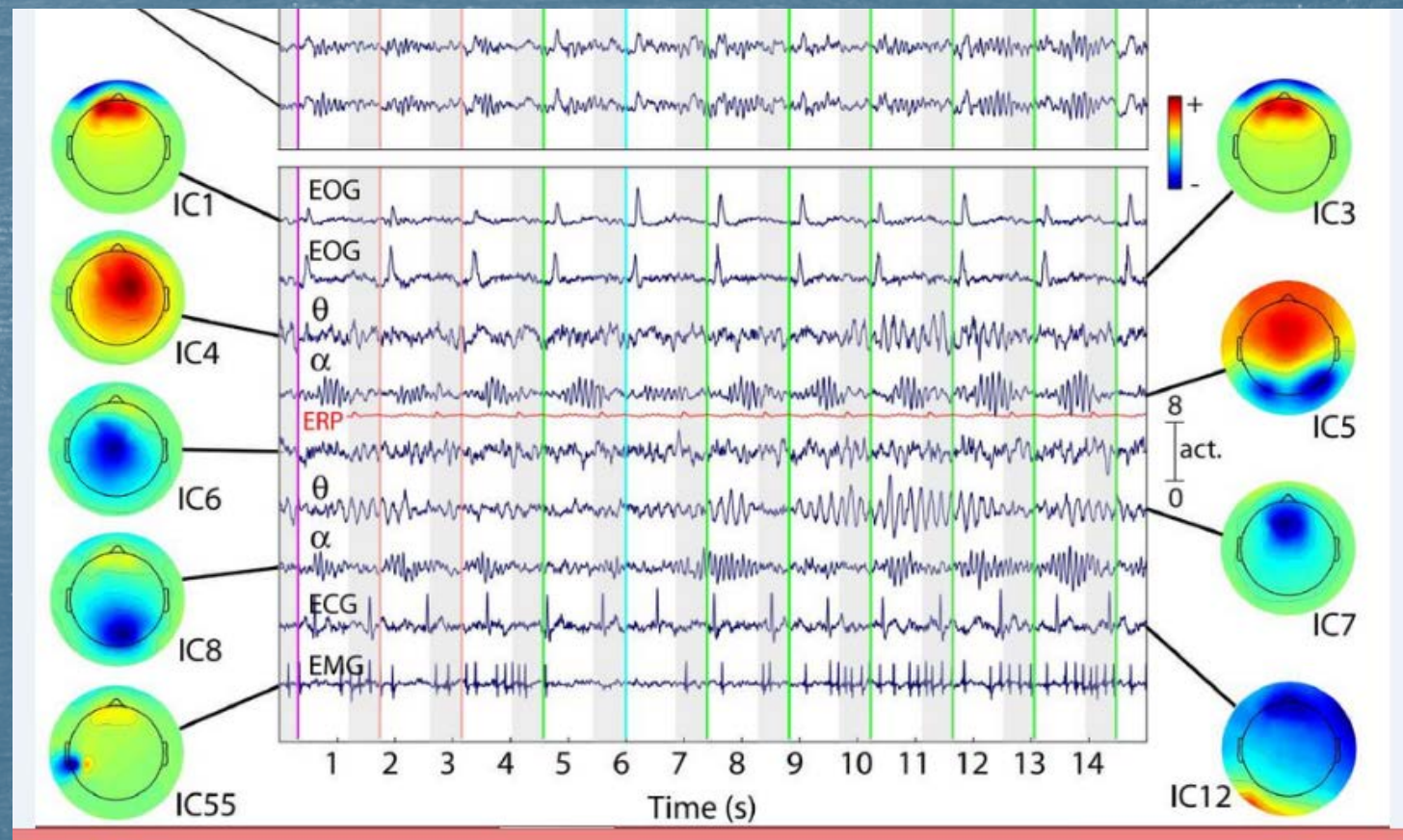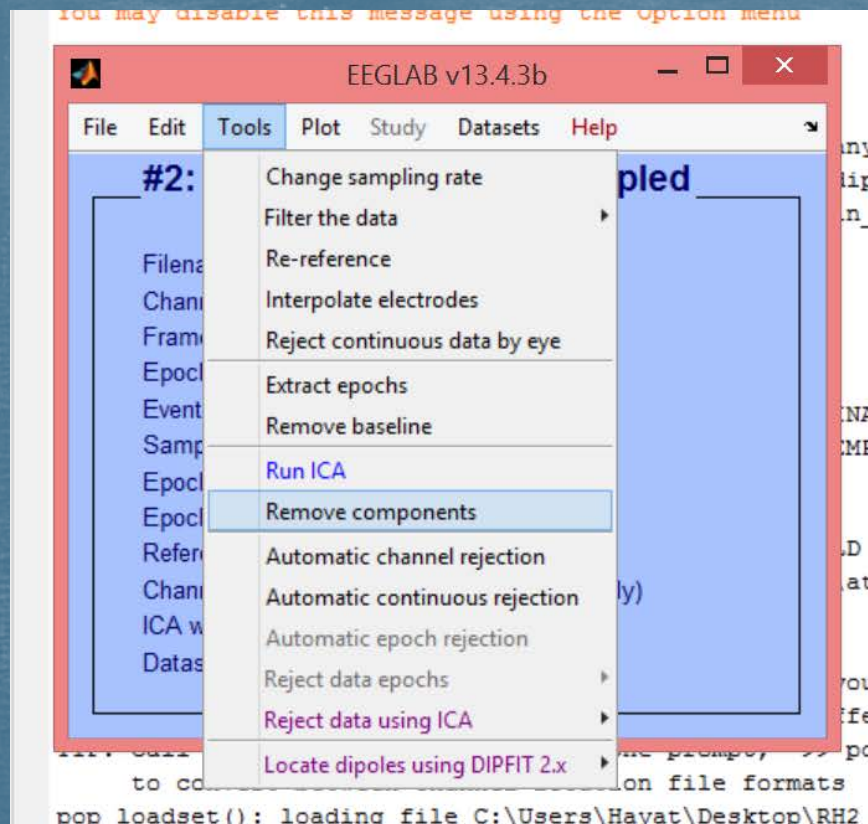

Removing artifacts (remove EOG, EMG and ECG)